

Chester Chest[™] – 2400/2402 and 2410/2412 USER MANUAL







Thank You For Your Purchase!

Thank you for your purchase of VATA's Chester Chest™.

Chester Chest™, an industry standard since 1987, enables physicians, nurses, home health staff, patients and support persons to develop competence with the most common types of vascular access devices. Within one realistic and portable training aid, is a great tool for teaching, training, competency testing and skills assessment. This is the most realistic and complete model for central line care offered anywhere! The area around all catheter sites can be used to practice cleansing, application of dressings and securement devices, fluid infusion and withdrawal.

Please read the entire User's Manual carefully before using the model to insure that you understand the proper care and use. This will also avoid situations that may not be covered by the warranty and help you enjoy the maximum benefit.

Contact Information

Whatever your question, problem, or comment, *VATA*'s Customer Service is here to help. There are four quick, easy ways to contact us, so you can choose what works best for you. If there's anything we can do, just let us know!

VATA Inc.

308 South Sequoia Parkway, Canby, Oregon 97013, USA

Phone: **503.651-5050**Fax: **503.651-5052**E-Mail: *info@vatainc.com*Web: *www.vatainc.com*



Overview of Chester Chest™

The left chest area of Chester Chest™ features:

- Outer Tissue Flap 0405/0404 A specially formulated material duplicates the feel of human tissue which, when placed over the chest port, provides a realistic practice of palpating and accessing, with proper access of the port being confirmed by a blood withdrawal – just like the real thing!
- · A real implanted port under the Outer Tissue Flap
- Difficult Accessing Inserts 0420/0421, 0430/0431 and 0440/0441 which are made
 of a soft tissue like material and placed either under or over the port to simulate
 palpating and accessing a port with one of the following types of placements:
 normal, "tipping", "wandering" or "deeply placed".
- A rigid underlying surface with molded ribs and a recessed area for the interchangeable inserts

The right chest area of Chester Chest™ features:

- A 9.6FR tunneled central catheter that is visible up to the clavicle. The Dacron cuff is also discernable.
- An external jugular vein is slightly raised with an opening for you to insert your own catheter and there is also an opening in the clavical area for the placement of a subclavian catheter.

The arm of Chester Chest™ features:

- A dual lumen 5FR PICC (Peripherally Inserted Central Catheter) exiting the basilic vein from the inner bicep area, which is currently the preferred site for placement.
 The basilic vein is slightly raised for easy identification.
- A pre-positioned 20g IV catheter in the forearm.

*Note: There is also an option to purchase Chester Chest with a Port Access Arm (pictured on page 7), which includes features of the standard arm plus a recessed area for the placement of a peripheral arm port. (Port not included – see part 0417 in Models & Parts on page 15). The base of the recessed area is made of soft material that permits the port to "float" when accessed. This area is covered with a removable tissue-like flap which, when placed over the port, provides the realistic feel of palpating and accessing.

All central lines are attached to a simulated blood reservoir bag in either the torso or the arm to permit the practice of "blood" withdrawal and fluid infusion. Cleansing, application of dressings and securement devices can be demonstrated at all sites.

3




Overview of *Chester Chest™*, cont.

Chester Chest™ can be used in either the upright or supine position. Chester is available in both lightly or darkly pigmented versions. If you would prefer, you can send us the catheters your institution uses and we can install them for you.

A video overview of Chester can be found on our website:

https://vatainc.com/resources/videos

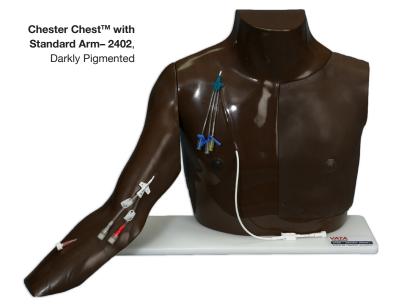



Chester Chest™ - 2400




- 1 0406 Real Port (IVAD)
- 2 0407 Tunneled Cental Catheter 9.6FR
- 3 Base
- 4 0405 Outer Tissue Flap
- **5 0451** Simulated Blood Reservior Bag for Torso
- 6 0430 Difficult Accessing Insert simulates "tipping" port
- **0420** Difficult Access Insert simulates a "wandering" port
- 8 0440 Difficult Access Insert simulates a deeply placed port
- 9 0458 Talc-Cornstarch

- IV Catheter 20G
- User's Manual
- 12 0408 Dual Lumen PICC 5FR
- 13 2388 Standard Arm
- **14 0410 Optional** Triple Lumen Catheter Shown (See Models & Parts on p.15)
- Opening for Optional Subcalvian Catheter
- Opening for Optional Jugular Catheter
- Torso




Backside of Chester Chest™ Torso

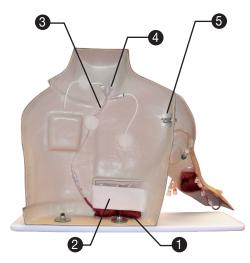

- 1 0451 Simulated Blood Reservoir Bag for Torso
- 2 Retention Strap for Simulated Blood Reservoir Bag
- **3 0450** Three Way Parallel Tubing Set for Torso
- 4 Extra connection for use with with Optional Triple Lumen, Subclavian and Jugular Catheters
- **5 0455** Bolt and Wing Nut to Attach Arm to Torso


Standard Arm

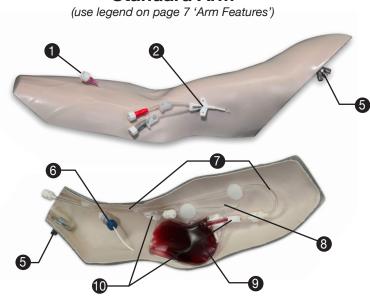



Port Access Arm

The Port Access Arm has the option to include an Optional Peripheral Port 0417 posterior to the PICC insertion site in the recessed area of the arm.

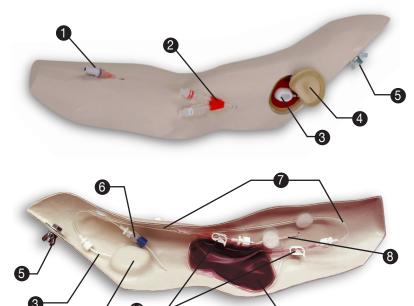

Arm Features

- 1 IV Catheter 20G
- 2 0408 Dual Lumen PICC 5FR
- **3 0417** Optional Peripheral Port (See Models & Parts on p.15)
- 4 0442 skin flap covering recessed area for optional peripheral port 0417
- **5 0455** Bolt and Wing Nut to Attach Arm to Torso

- 6 Dual PICC Attachment to Tubing Set
- **0446** Tubing Set for Arm
- 8 20G IV Catheter Attachment to Tubing Set
- 9 0453 Simulated Blood Reservoir Bag for Arm
- 10 Snap Clamps
- Back of Recessed Area

https://vatainc.com/resources/videos

Setting up your Chester Chest™



Unwrap the arm and body. Set the bag of additional parts, packed with the torso, aside for use later. Place the torso on a flat surface in an upright position. Remove the wing nut from the arm, keeping the bolt in place, and align the bolt with the hole in the right shoulder; insert and reattach the wing nut (see fig. 1).




If you are going to use injection caps or needleless access caps, attach before attaching the simulated blood reservoir bag. On the backside of the torso, locate the simulated blood reservoir bag and remove from torso. Carefully remove the white cap on the simulated blood reservoir bag, taking care not to let the liquid leak out and attach the female leur fitting to the male leur fitting on the bottom of the three-way parallel tubing (see *fig. 2*).



Replace the reservoir bag behind the white retention strap, making sure that the tubing attached to the bag is positioned exiting the bag from the bottom (see *fig.* 3).




Do not position reservoir bag with tubing exiting from the top as you will withdraw air (see fig. 4).

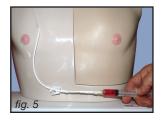

To prime the tubing with simulated blood attach a syringe to the tunneled central venous catheter and drawback to remove all the air in the line. Be sure the snap clamp is open on the catheter. This procedure may need to be repeated depending on the size of your syringe. You are done when the simulated blood is visible in the syringe (see fig. 5).






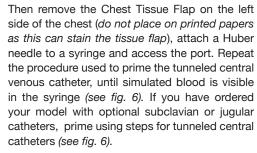




The torso is now primed and ready for use. Please note that one of the three parallel tubing sets is not used and has a cap attached to the male leur connector. This tubing set is provided should you wish to add a subclavian or jugular catheter (purchase of a special connector may be needed for some catheters).

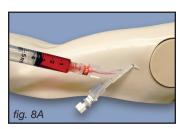

To prime the arm tubing, attach injection caps or needless access devices on the IV catheters and dual lumen PICC before attaching the simulated blood reservoir bag. On the backside of the arm, open all snap clamps (see fig. 7).



Attach a syringe to each side of the dual lumen PICC and withdraw air, until simulated blood is visible in the syringe. This procedure may need to be repeated, depending on the size of your syringe. Repeat this procedure with the IV catheter (see fig. 8A & 8B). If you have ordered your arm with the **optional** peripheral port, attach a Huber needle to a syringe, access the port and withdraw air until simulated blood is visible in the syringe (see fig. 9). The arm is now primed and ready for use.

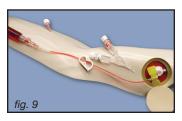



Proper Use of your Chester Chest™

Chest Tissue Flap - 0405

When removing the Chest Tissue Flap, always pull gently from the edge to prevent damage to the flap. Due to the elasticity of the Chest Tissue Flap, it may be easier to attach when Chester Chest™ is in a supine position. In order to ensure a realistic feel when palpating the port (IVAD), the Chest Tissue Flap is formulated to be soft to the touch. As the material is soft, care must be taken to maximize the useful life. There are three areas that should be avoided, as all can cause premature tearing of the Chest Tissue Flap:



Do not pull the Chest Tissue Flap back to view the placement of the needle in the port septum (see fig. 10).




Do not press fingernails into the tissue flap when palpating the port (see fig. 11).





Do not "rock" the access needle back and forth to confirm placement (see fig. 12).



fig. 13B

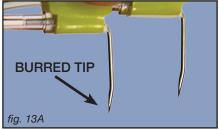

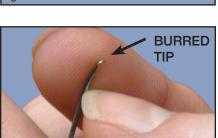

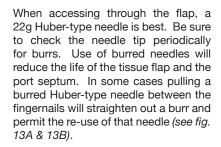

Should you need additional Huber practice needles you can order:

5025 3/4" 20G right angle Huber needle 8" extension set 1 dozen.

5026 1" 20G right angle Huber needle 8" extension set 1 dozen.



The entire area of the tissue flap can be used for accessing. Just move the port to the desired location and place the Chest Tissue Flap over it (see fig. 14). See text under **Difficult Accessing Insert 0430** to learn how to extend catheter length to move port (see fig. 18).



Cleaning of the Chest Tissue Flap or Difficult Accessing Inserts is best done by wiping with an alcoholwetted, non-linting cloth. The part should be permitted to air dry and then powdered with the cornstarch talc supplied. Any excess talc can be removed with a dry cloth. If at any time the skin flap becomes tacky to the touch, talc should be applied (see fig. 15).



Difficult Accessing Inserts


The model is supplied with three different Difficult Accessing Inserts. These are used alone or in conjunction to simulate various difficult accessing scenarios. All the inserts have the part numbers on them for easy identification.



0420 Difficult accessing insert simulates a "wandering" port



0430 Difficult accessing insert simulates "tipping" port

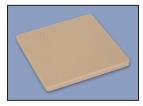

0440 Difficult accessing insert simulates a deeply placed port





DIFFICULT ACCESSING INSERT 0420 for simulation of a "wandering" or "floating" port. Remove the 0430 insert and replace with the 0420. The #420 should be located in the upper left hand corner. For best results in simulating a wandering or floating port, place a small amount of K-Y or other lubricating jelly under and on top of your port. This will vary by the type of port used (see fig. 16). Place the port in the center of the depressed area and reattach the Chest Tissue Flap. When you palpate the port, it will move around.







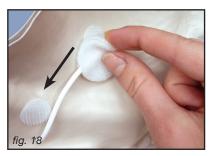

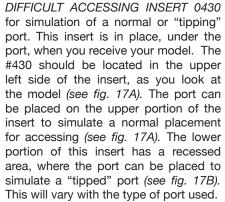

If you find you need additional catheter length to place the port in this area, look on the backside of the torso where the port catheter comes through from the front side. There is Velcro on the catheter tubing to secure the catheter to the torso. By pulling the Velcro off the upper placement and moving it to the lower attachment position, you will get an additional length of tubing (see fig. 18).



DIFFICULT ACCESSING INSERT 0440 for simulation of a "deeply placed" port. The #440 will be located in one of the corners. This piece is placed over the top of the port, with the Chest Tissue Flap then placed over the 0440. Best results using the 0440 insert are experienced when the 0420 insert is used in the recessed area under the port (see fig. 19).



Chester Chest™ - Refurbishing Service


The purchase of training models is a significant financial investment. Realizing this, **VATA** is proud to offer a refurbishing service to bring your well-used **Chester Chest** back to life for about half the cost of a new model!

How It Works:

14 -

- Contact VATA to receive your Return of Goods Authorization (RGA) number and ship your Chester to our facility.
- Upon receipt, your Chester Chest will be evaluated for needed parts and repairs. We will create an estimate for your approval and prepayment is required prior to starting the refurbishment process.
- The maximum refurbishing cost is \$450 + shipping* (as long as you have the body shell, arm shell, and base).
- Refurbishment takes two to five business days to complete after we have received your approval and prepayment for the work to be completed.
- Return shipping is typically around 5 business days (the cost of shipping will be included in the estimate).

To begin the process, please call 503-651-5050 or email info@vatainc.com to receive your Return of Goods Authorization (RGA) number.

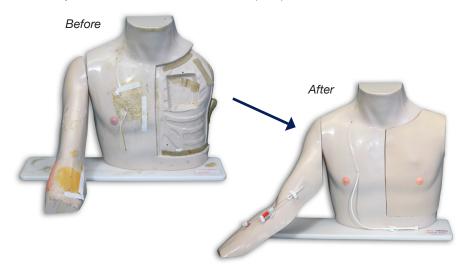

*This refurbishment pricing is valid January 1, 2021-December 31, 2021.



Chester Chest™ - Models & Parts

The following items are available for purchase at www.vatainc.com.

All parts on this model are available individually.

Chester Chest™

2400	Chester Chest™ with Standard Arm, lightly pigmented
2402	Chester Chest™ with Standard Arm, darkly pigmented
2410	Chester Chest™ with Port Access Arm, lightly pigmented
2412	Chester Chest™ with Port Access Arm, darkly pigmented

Accessories & Parts

0401	Optional Carrying Case for Chester Chest™
0404	Chest Tissue Flap Replacement, darkly pigmented for 2402/2412
0405	Chest Tissue Flap Replacement, lightly pigmented for 2400/2410
0406	Practice Port (IVAD)
0407	Tunneled Central Venous Catheter 9FR (CVC)
0408	Dual PICC Catheter 5 FR (requires part 0418 for attachment)
0409	Optional Tunneled Dual Lumen Catheter 9.6FR
0410	Optional Triple Lumen Catheter (requires part 0454 for attachment)
0417	Optional Peripheral Arm Port for Port Access Arm
0418	6 FR Universal Catheter Connector
0442	Replacement Skin Flap to cover peripheral arm port, lightly pigmented
0443	Replacement Skin Flap to cover peripheral arm port, darkly pigmented
0446	Clear Tubing for Chester Arm
0450	Three-way Parallel Tubing Set
0451	Simulated Blood Reservoir Bag for Chester torso
0453	Simulated Blood Reservoir Bag for Chester arm
0454	9.6 FR Universal Catheter Connector
2386	Replacement Port Access Arm for Chester, lightly pigmented
2387	Replacement Port Access Arm for Chester, darkly pigmented
2388	Replacement Standard Arm for Chester, lightly pigmented
2389	Replacement Standard Arm for Chester, darkly pigmented
1491	Simulated Blood, One Quart - New, Stain Resistant
1494	Simulated Blood, One Gallon - New, Stain Resistant
5025	3/4" 20G right angle Huber needle (set of 1 dozen)
5026	1" 20G right angle Huber needle (set of 1 dozen)



You Might Also Like These Products...



Vascular Access Ultrasound Phantom, Lightly Pigmented – 0705 Vascular Access Ultrasound Phantom, Darkly Pigmented – 0706



Port – Body in a Box™, Lightly Pigmented – 5010 Port – Body in a Box™, Darkly Pigmented – 5011





Advanced Venipuncture Training Aid™, Lightly Pigmented – 2365
Advanced Venipuncture Training Aid™, Darkly Pigmented – 2366

